PROMESH® SURG ABSO VENT

SEMI RESORBABLE PARIETAL REINFORCEMENT IMPLANT EXTRAPERITONEAL

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CE 0297

NOT304_260929
Version du
29/09/2026

℞ only Caution : Federal law (USA) restricts this device to sale, distribution and use by or on the order of a physician.

This release is the last update of the instructions of use and replace the previous edition

Product under regulation 2017/745 (EU).

PROMESH® SURG ABSO VENT

SEMI RESORBABLE PARIETAL REINFORCEMENT IMPLANT
EXTRAPERITONEAL
SINGLE USE STERILE PRODUCT

MANUFACTURER IDENTIFICATION

COUSIN BIOTECH
ALLEE DES ROSES
59117 WERVICQ SUD
FRANCE

IDENTIFICATION AND DESCRIPTION OF THE DEVICE

PROMESH® SURG ABSO VENT meshes are extraperitoneal semi-resorbable parietal reinforcement implants.

INTENDED USE

Intended to be used as parietal reinforcement implant.

INDICATIONS

Repair of ventral hernia that require the addition of an extraperitoneal reinforcing or bridging material to obtain the desired surgical result.

CONTRAINDICATIONS

Do not use in the following cases:

- Allergy to any of the components
- Infected site
- Pregnancy
- Growing children

UNDESIRABLE SIDE EFFECTS

Like any implantable medical device, this implant is susceptible to generate possible undesirable side effects that may lead to reoperation:

- Discomfort/Pain
- Infection
- Inflammation
- Fistula
- Chronic Pain
- Recurrence
- Adhesion
- Erosion/Extrusion
- Seroma/lymphocele
- Hematoma
- Mesh migration
- Mesh deformation
- Allergic reaction
- Foreign body reaction
- Irritation nearby organ
- Ileus / Obstruction
- Organ injury

TARGET POPULATION

Adult male or female suffering from a symptomatic ventral hernia necessitating a surgical repair with a method that includes an extraperitoneal implantation.

IMPLANTED MATERIALS

Polypropylene (PP) (non-resorbable) – Poly L Lactic Acid (PLLA) (resorbable)

	Weight PP (g)	Weight PLLA (g)
PSABVT12DI	0,70	1,05
PSABVT1515	1,39	2,08
PSABVT1530	2,78	4,18
PSABVT2025	3,09	4,64
PSABVT3030	5,57	8,36

Origin neither human nor animal – Semi resorbable.

PERFORMANCES

The device PROMESH® SURG ABSO VENT contributes to reduce symptoms from hernia and contributes to avoid complications of hernia.

LINK TO SUMMARY OF SAFETY AND PERFORMANCE CHARACTERISTICS:

https://www.cousin-biotech.com/d-6318526d8b749_under-construction.pdf

MRI SAFETY

The implant has not been tested for safety and compatibility in the case of postoperative investigation using Magnetic Resonance Imaging (MRI). However, the implant is mainly made with Polypropylene and Poly-L-Lactic Acid: non-metallic components. Consequently, PROMESH® SURG ABSO VENT is considered MR safe

PRECAUTIONS FOR USE

PROMESH® SURG ABSO VENT must always be separated from abdominal cavity by peritoneum
PROMESH® SURG ABSO VENT implants are delivered sterile (ethylene oxide sterilisation).
When opening the carton box, the surgeon has to make sure that the lot number, the model and the size of the mesh are equal to the data on the label of the inner pouch. Otherwise, the use of the prosthesis is strictly prohibited.
Before any use, inspect the integrity of packaging and device (of which peelable pouches). Do not use in the event of deterioration of the device and/or the packaging. Do not use if the device is out of date. COUSIN BIOTECH does not offer any guarantee or recommendation as far as the use of a particular type of means of fixation is concerned. Mesh fixation points should be at least 1cm from the edge of the mesh with 1cm spacing between fixation points.
Unpacking and manipulation of the implant must be carried out following asepsis standards in order to guarantee the sterility.
Prevent the mesh from any contact with objects that may alter their surface.
This device must be implanted by a qualified surgeon (familiar with the relevant anatomy and experienced in visceral surgery).

IMPORTANT: DO NOT REUSE - DO NOT RESTERILIZE

As specified on the product's labelling, the **PROMESH® SURG ABSO VENT** mesh is for single use only. It can not be re-used and/or re-sterilized (potential risks would be and are not limited to: loss of the product's sterility, risk of infection, loss of the product's efficiency, recurrence)

INSTRUCTIONS TO BE GIVEN TO THE PATIENTS BY THE SURGEON

The surgeon should inform the patient of the potential physical and psychological restrictions and consequences of implanting the device. The patient must be informed of the surgical risks and possible side effects. The surgeon should invite the patient to return for further consultation if the patient presents symptoms that appear abnormal. An implant card for the patient is provided by the healthcare professional at the clinic. This card provides information to identify the device and traceability elements as well as the name, address and website of the manufacturer. Please refer to the section “Implant card information leaflet” for instructions on how to complete it. **The surgeon shall invite the patient to scan the implant card immediately upon receipt to keep track of it in case of loss**. The summary of the safety and performance characteristics of the device and the instructions for use can also be found on the COUSIN BIOTECH website.

STORAGE PRECAUTIONS


To be stored in a dry place away from sunlight and at room temperature in its original packaging.

SURGICAL TECHNIQUE

Extra-Peritoneal Hernia repair

Example: Retrorectus mesh repair of a ventral hernia (Rives-Stoppa technique).
This technique utilizes the hernia sac to separate the mesh from the intra-abdominal contents.

- Superior to the umbilicus, dissection is performed above the posterior rectus fascia and under the rectus muscle.
- Below the umbilicus, dissection occurs in the preperitoneal space due to the lack of a posterior rectus sheath.
- The PROMESH® SURG ABSO VENT implant is placed in the newly formed pre-peritoneal space. The PROMESH® SURG ABSO VENT mesh has an oriented elasticity:

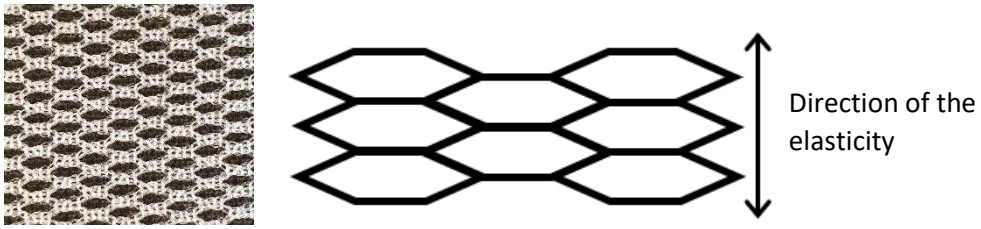- Fixation of the mesh if needed, to the muscle layer above
- The peritoneum must be meticulously checked and no defect in the peritoneum should be left open to prevent any contact between the mesh and intra-abdominal contents.
- Redon’s drains under aspiration are usually placed. The anterior sheath is closed. If tension is present, relaxing incisions may be used
- Skin closure.
- Using the Redon’s drains, fill the pre-peritoneal space with an antiseptic solution and clamp the drains for 1 hour.
- A belt can be used to wrap the patient’s abdomen for a few days to reinforce the abdominal wall.

EXPLANTATION AND DISPOSAL OF DEVICES

This device is intended for permanent parietal reinforcement and to remain permanently implanted. Routine removal or replacement is not required. In the event of certain complications, a secondary surgical procedure for partial or complete mesh explantation may be considered, if clinically indicated. Complete removal of all implant material may not always be technically feasible. The decision to explant the device lies with the responsible surgeon and shall be based on the individual clinical situation and current medical practice.
Mesh explantation is a technically demanding reoperation and should only be performed by a surgeon who is appropriately qualified, trained, and experienced in mesh removal techniques. As with any reoperation, mesh explantation carries significant risks, including but not limited to bleeding, infection, injury to adjacent anatomical structures (e.g. nerves, vessels, bowel, bladder, or other organs), anaesthesia related complications, recurrence or persistence of the original condition, worsening of symptoms, chronic pain, and the need for additional surgical interventions. Patients should be informed that mesh explantation may resolve, improve, partially improve, fail to resolve or worsen the reported symptoms or adverse effects.
If a hernia mesh implant is explanted and retrieval data are to be collected, the device should be retrieved and handled in accordance with ISO 12891 1:2015 “Implants for surgery – Retrieval and analysis of surgical implants – Part 1: Retrieval and handling”.

Explanted devices that are not returned for analysis must be disposed of in accordance with applicable local regulations for infectious and hazardous clinical waste.

IMPORTANT

For more information about the use of this product, please contact your COUSIN BIOTECH representative or distributor.

INFORMATION REQUEST AND COMPLAINTS

Following its quality policy, COUSIN BIOTECH is committed to making every effort to produce and supply a high-quality medical device. However, if a health professional (client, user, prescriber, etc.) has a complaint or cause for dissatisfaction with a product in terms of quality, safety or performance, they must inform COUSIN BIOTECH as soon as possible. In case of failure of an implant or if it has contributed to causing a serious adverse reaction in the patient, the health centre must follow the legal procedures in that country and inform COUSIN BIOTECH immediately. For any correspondence, please specify the reference, batch number, contact person's details and a comprehensive description of the incident or complaint.
Brochures, documentation and surgical technique are available on request from COUSIN BIOTECH and its distributors. If further information is needed or required, please contact your COUSIN BIOTECH representative or distributor.

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.


IMPLANT CARD INFORMATION LEAFLET

How to complete the implant card ?

The implant card, for the patient, should be completed as follows:

Front side of the implant card

The user has to complete the implant card before transmitting the implant card to the patient. Please complete it as follows:

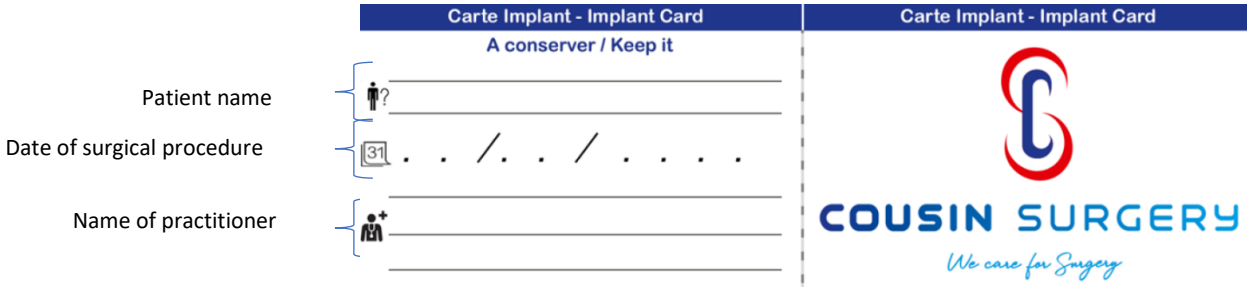Back side of the implant card

No fields need to be completed by the user. The following information on the device is printed directly by the manufacturer on the back of the implant card. A QR code allowing the patient to consult the instruction for use, as well as the website address, are available on the implant card.

Information to identify the device	Traceability information
• Indication that the product is a medical device • Description of the device family • Trade name of the device • Reference of the device • Device is MR safe • Denomination name of the device in sales countries	• UDI number • Batch number • Name and address of the legal manufacturer

Explanation of the used symbols in the Implant Card

All symbols are commonly used and explained at the end of the instruction for use.

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